PERSONALIZED COMFORT DRIVER SEAT WITH INFLATABLE LUMBAR PAD, HEAT & MASSAGE

FOR USE AT HOME, IN THE OFFICE OR IN YOUR VEHICLE

Please read this manual thoroughly before operating your seat cushion. This manual should be saved for future reference.

Materials & Electrical Components

- 100% polyester exterior; polyurethane foam interior
- CA 117 fire retardant
- Car / cigarette lighter adapter
- AC adapter
- Hand-held remote

OPERATING YOUR OBUSFORME SEAT CUSHION

Using the Car Adapter

Plug the car adapter into your vehicle's cigarette lighter or auxiliary power source. The car adapter requires 12 volts.

Warning! The car adapter must NOT be plugged into any cigarette lighter or auxiliary power source until you have confirmed that the cigarette lighter or auxiliary power source does not exceed 12 volts. Plugging the car adapter into the cigarette lighter or auxiliary power source found in some vehicles may cause the seat cushion to overheat and malfunction.

The car adapter contains a three-amp fuse (F3A). If the car adapter is used properly, this fuse should not burn out. If it does, it can be replaced with another F3A fuse which can be purchased at most hardware or electronics stores. To replace the fuse, simply unscrew the top portion of the car adapter. Remove the old fuse and replace it with a new F3A.

Warning! Replace the F3A in the car adapter ONLY with another F3A. Using any other fuse may cause the car adapter to overheat and malfunction. If you are uncertain as to which fuse to use, take the old fuse to your local hardware or electronics store and consult a professional. Plugging in the car adapter will automatically activate the heat function. To turn the heat off, simply unplug the car adapter (or use the hand-held remote.)

Warning! When not in use, the seat cushion car adapter should be unplugged from the cigarette lighter or auxiliary power source. If the car adapter is not unplugged and you remain seated on the cushion, it will continue to drain the battery of your vehicle even if your vehicle is not running.

Using the AC Adapter

Connect the AC adapter plug into a power source of the same voltage, frequency and rating as indicated on the transformer. Activate and control the desired functions using the hand-held remote.

How to use the Massage Option

The massage and pulse options can be applied to your upper back, lower back or thighs, either together or independently. To activate the massage, adjust the "UPPER BACK", "LOWER BACK" and "THIGHS" switches to your preferred settings. Control the intensity of your massage/pulse with the "LOW / MEDIUM / HIGH" switch.

How to use the Heat Option

To activate the heat using the hand-held remote*, push the "HEAT" button once. The heat indicator light will illuminate. To turn the heat off, depress the "HEAT" button a second time. The heat indicator light will go out. Please see "Using the Car Adapter" for important safety information.

SAFETY DEVICES

For your protection, these seat cushions have three automatic safety devices:

1 30-Minute Automatic Shut-off - Massage

After 30 minutes, the massage and pulse functions will automatically shut off. You can turn these features on again after they have shut off. The heat function DOES NOT automatically shut off after 30 minutes. To turn off the heat function, you must press the "HEAT" button or unplug the car adapter from the power source.

2 Weight-Sensitive Automatic Shut-off

The ObusForme Seat Cushions are built with a pressure switch that will automatically turn off the heat function when your body weight is removed from your ObusForme Seat Cushion.

3 Temperature-Regulated Automatic Shut-off

The ObusForme Seat Cushions are designed with an internal thermostat that regulates the temperature of the cushion. Your cushion will automatically shut off when the temperature reaches its maximum heat level and be reactivated once the temperature lowers.

Limited 1 Year Warranty

HoMedics sells its products with the intent that they are free of defects in manufacture and workmanship for a period of one year from the date of original purchase, except as noted below. HoMedics warrants that its products will be free of defects in material and workmanship under normal use and service. This warranty extends only to consumers and does not extend to Retailers. To obtain warranty service on your HoMedics product, mail the product and your dated sales receipt (as proof of purchase), postpaid, to the following address:

HoMedics Group Canada
ATTN : Consumer Relations,
344 Consumers Road
Toronto, Ontario M2J 1P8 CANADA
No COD's will be accepted.

Phone: (416) 785-1386
Fax : (416) 785-5862
Toll Free 1-888-225-7378
8:30 a.m. to 5:00 p.m. ET
Monday – Friday

HoMedics does not authorize anyone, including, but not limited to, Retailers, the subsequent consumer purchaser of the product from a Retailer or remote purchasers, to obligate HoMedics in any way beyond the terms set forth herein. This warranty does not cover damage caused by misuse or abuse; accident; the attachment of any unauthorized accessory; alteration to the product; improper installation; unauthorized repairs or modifications; improper use of electrical/power supply; loss of power; dropped product; malfunction or damage of an operating part from failure to provide manufacturer's recommended maintenance; transportation damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of HoMedics. This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adoption to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.


THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. THERE SHALL BE NO OTHER WARRANTIES EXPRESS OR IMPLIED INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS OR ANY OTHER OBLIGATION ON THE PART OF THE COMPANY WITH RESPECT TO PRODUCTS COVERED BY THIS WARRANTY. HOMEDICS SHALL HAVE NO LIABILITY FOR ANY INCIDENTAL, CONSEQUENTIAL OR SPECIAL DAMAGES. IN NO EVENT SHALL THIS WARRANTY REQUIRE MORE THAN THE REPAIR OR REPLACEMENT OF ANY PART OR PARTS WHICH ARE FOUND TO BE DEFECTIVE WITHIN THE EFFECTIVE PERIOD OF THE WARRANTY. NO REFUNDS WILL BE GIVEN. IF REPLACEMENT PARTS FOR DEFECTIVE MATERIALS ARE NOT AVAILABLE, HOMEDICS RESERVES THE RIGHT TO MAKE PRODUCT SUBSTITUTIONS IN LIEU OF REPAIR OR REPLACEMENT.

This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or sales of such products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate as to any products or parts thereof which are repaired, replaced, altered, or modified, without the prior express and written consent of HoMedics. This warranty provides you with specific legal rights. You may have additional rights which may vary from country to country. Because of individual country regulations, some of the above limitations and exclusions may not apply to you.

For more information regarding our product line in Canada, please visit: www.homedics.ca

HOW TO SECURE YOUR OBUSFORME SEAT CUSHION TO A VEHICLE SEAT

1. Slide strap "A" and strap "B" over the back of your vehicle seat. Adjust the seat cushion on your vehicle seat to your personal preference.
2. Attach strap "LC" to strap "RC", securing with the plastic hooks located on the end of each strap.
3. Attach strap "LD" to strap "RD", securing with the plastic hooks located on the end of each strap.
4. Pull strap "C" over your headrest and secure it to strap "A" using the buckle attached.

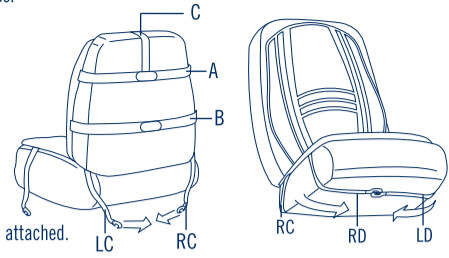Warning!

Use caution when stretching and securing the elastic straps. The elastic straps may recoil and cause injury if you fail to properly secure the plastic hooks before releasing the elastic straps.

CARING FOR YOUR OBUSFORME SEAT CUSHION

All servicing of your seat cushion MUST be performed by authorized ObusForme personnel. Please see warranty information for further details, or contact ObusForme for instructions on returning it for examination.

To Clean

Wipe the fabric cover with a damp cloth or sponge and a mild detergent. To dry, pat with a clean, dry cloth. DO NOT rub the seat cushion excessively, immerse in any liquid or machine wash.

DO NOT use abrasive cleaners, gasoline, kerosene, glass/furniture polish, paint thinner, or water above 49°C (120°F) to clean your seat cushion.

Warning!

Before cleaning your seat cushion, ensure that the car adapter or AC adapter is unplugged from the power source. Do not plug the car adapter or AC adapter into any power source until your seat cushion is completely dry.

To Store

Fold and return your ObusForme Seat Cushion to its box and/or store in a cool, dry place.

OTHER IMPORTANT SAFETY TIPS

Your seat cushion is intended for indoor* or vehicle use ONLY.

CAUTION! To reduce the risk of burns, fire, electric shock or injury:

- Always unplug your ObusForme Seat Cushion from the cigarette lighter, auxiliary power source or electrical outlet immediately after using and before cleaning.
- Do not place or store your ObusForme Seat Cushion where it can fall into water or any other liquid.
- Only use your ObusForme Seat Cushion for its intended use, as described in this manual.
- Do not use any attachment not recommended by ObusForme.
- Never operate your ObusForme Seat Cushion if the cord or plug is damaged, if it is not working properly or if it has been damaged in any manner.
- Use the car adapter or AC adapter provided by ObusForme. Do not substitute any other car adapter or AC adapter; doing so will void the product warranty.
- Disconnect power when not in use.
- Do not use any cord as a handle.
- Keep all cords away from heated surfaces.
- Close supervision is necessary when children or persons with disabilities use the ObusForme Seat Cushion.
- Never leave the ObusForme Seat Cushion unattended when plugged in.

CAUTION! Consult your doctor before using an ObusForme Seat Cushion if you have any concerns regarding your health, including the following:

- if you have a pacemaker
- if you are experiencing unexplained pain or swelling
- if you've recently experienced serious muscular injury
- if you are, or may be, pregnant
- if you have sensitive skin

Do not use your seat cushion as a substitute for medical attention.

Discontinue use of your seat cushion if you experience any skin inflammation or skin eruptions.

ASSOCIATION OF CANADIAN ERGONOMISTS

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