BACK MASSAGER WITH HEAT

INSTRUCTION MANUAL AND WARRANTY INFORMATION

ONE YEAR LIMITED WARRANTY

BB-BMS

IMPORTANT SAFETY INSTRUCTIONS

WHEN USING ELECTRICAL PRODUCTS, ESPECIALLY WHEN CHILDREN ARE PRESENT, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED, INCLUDING THE FOLLOWING:

READ ALL INSTRUCTIONS BEFORE USING

DANGER – TO REDUCE THE RISK OF ELECTRIC SHOCK:

- Always unplug the appliance from the electrical outlet immediately after the massager is fully charged.
- DO NOT reach for an appliance that has fallen into water. Unplug it immediately.
- DO NOT use while bathing or in a shower.
- DO NOT place or store appliance where it can fall or be pulled into a tub or sink. DO NOT place in or drop into water or other liquid.
- Never use pins or other metallic fasteners with this appliance.

WARNING – TO REDUCE THE RISK OF BURNS, FIRE, ELECTRIC SHOCK, OR INJURY TO PERSONS:


- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by HoMedics; specifically any attachments not provided with the unit.
- NEVER operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to HoMedics Service Center for examination and repair.
- Keep the cord away from heated surfaces.
- DO NOT operate where aerosol (spray) products are being used or where oxygen is being administered.
- DO NOT operate under a blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- DO NOT use outdoors.
- Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
- Do not use this massager in close proximity to loose clothing, jewelry, or long hair.
- Never operate on a soft surface such as a bed or couch where the air openings may be blocked.
- This product has an internal, non-replaceable lithium battery. This battery is not user serviceable. The battery must be removed from the appliance before it is scrapped and the appliance must be disconnected from the supply mains when removing the battery.

Please dispose in accordance with local, state, province, and country regulations.


WARNING: For the purposes of recharging the battery, only use the detachable supply unit provided with this appliance.

BACK MASSAGER WITH HEAT

1. Your massager may arrive partially charged or with no charge. It is recommended that the user fully charges the massager before using. The massager is fully charged after approximately 4.5 hours.
2. Detach the Velcro panels and place back massager around the lower back. Adjust the Velcro to ensure the massager is snug around your back and the massage controller is in the front.

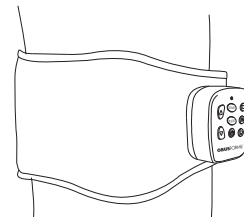3. Press the power button to turn the unit on. Press the vibration button to start the massaging function. Press the heat button to start the heating function, press again for higher heat. Press the Air compression button to start this feature. Note: all 3 functions can run at the same time. Each function needs to be turned on for it to work.

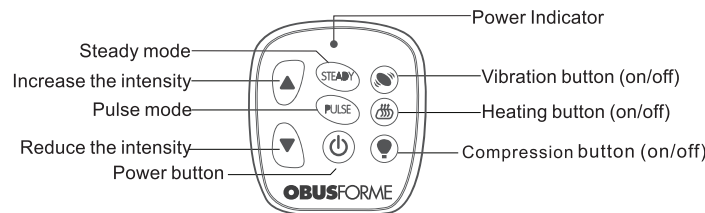Caution:

All servicing of this massager must be performed by authorized ObusForme service personnel only.

HoMedics Group Canada
6460 Kennedy Road, Unit C
Mississauga, Ontario,
Canada L5T 2X4

Telephone: 1-888-225-7378
8:30 a.m. to 5:00 p.m.,
Mon - Fri EST

www.ObusForme.com

LIMITED WARRANTY (FOR CANADA ONLY)

ObusForme guarantees this product will be free from defects in workmanship & materials for a period of one (1) year from the date of original purchase. This guarantee applies when the product is used for the purpose intended. In the advent of a defect as described above, please return this product together with the original purchase receipt to ObusForme for refund or replacement. This warranty gives you rights that vary by province/state. This warranty may change.

WHAT IS NOT COVERED

Wear & tear, aging (including foam discoloration, flattening, density, consistency), accidental damages, alterations, mishandling, faulty adjustment, misuse, improper care, power damage, rental use, discontinued items, service by anyone other than ObusForme™, shipping damages, neglect, items sold 'as is' or damage due to natural acts are NOT covered.

HOW TO OBTAIN WARRANTY SERVICE

Please obtain a Return Authorization Number & instructions prior to sending your item or it may be denied.

Please contact Customer Service by:

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