**OBUSFORME®**  
ENGINEERED FOR THE HUMAN BODY™

# SHIATSU & VIBRATION NECK MASSAGER WITH HEAT

INSTRUCTION MANUAL AND WARRANTY INFORMATION

**TWO YEAR LIMITED WARRANTY**

SM-SNM-02

## IMPORTANT SAFEGUARDS

WHEN USING ELECTRICAL PRODUCTS, ESPECIALLY WHEN CHILDREN ARE PRESENT, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED, INCLUDING THE FOLLOWING:

WHEN USING ELECTRICAL PRODUCTS, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED, INCLUDING THE FOLLOWING:

READ ALL INSTRUCTIONS BEFORE USING

**DANGER – TO REDUCE THE RISK OF ELECTRIC SHOCK:**

- Always unplug the appliance from the electrical outlet immediately after using and before cleaning.
- DO NOT reach for an appliance that has fallen into water. Unplug it immediately.

**WARNING – TO REDUCE THE RISK OF BURNS, ELECTRIC SHOCK, OR INJURY TO PERSONS:**

- An appliance should NEVER be left unattended when plugged in. Unplug from outlet when not in use and before putting on or taking off parts or attachments.
- Close supervision is necessary when this appliance is used by, on or near children, invalids or disabled persons.
- Use this appliance only for its intended use as described in this manual. DO NOT use attachments not recommended by ObusForme; specifically any attachments not provided with the unit.
- NEVER operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return it to ObusForme Service Center for examination and repair.
- Keep cord away from heated surfaces.
- NEVER drop or insert any object into any opening.
- DO NOT operate where aerosol (spray) products are being used or where oxygen is being administered.
- DO NOT operate under a blanket or pillow. Excessive heating can occur and cause fire, electrocution or injury to persons.
- DO NOT carry this appliance by the power cord or use cord as handle.
- To disconnect, turn all controls to the “off” position, then remove plug from outlet.
- DO NOT use outdoors.

- NEVER operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.
- NEVER operate on a soft surface such as a bed or couch where the air openings may be blocked.
- Use heated surfaces carefully. May cause serious burns. DO NOT use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.

## SAVE THESE INSTRUCTIONS

### CAUTION – PLEASE READ ALL INSTRUCTIONS CAREFULLY BEFORE OPERATING.

- This product is not intended for medical use. It is intended only to provide a luxurious massage.
- Consult your doctor prior to using this product, if
  - You are pregnant
  - You have a pacemaker
  - You have any concerns regarding your health
- Not recommended for use by Diabetics.
- Never leave the appliance unattended, especially if children are present.
- Never cover the appliance when it is in operation.
- Do not use this product for more than 15 minutes at a time.
- Extensive use could lead to the product's excessive heating and shorter life. Should this occur, discontinue use and allow the unit to cool before operating.
- Never use this product directly on swollen or inflamed areas or skin eruptions.
- DO NOT use this product as a substitute for medical attention.
- The product is equipped with overheating protection circuit, if the power is cut-off, please wait until the product is cooled down.
- DO NOT use this product before bed. The massage has a stimulating effect and can delay sleep.
- NEVER use this product while in bed.
- This product should NEVER be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls.
- This unit should not be used by children or invalids without adult supervision.
- NEVER use this product in automobiles.

#### **Caution:**

**All servicing of this massager must be performed by authorized ObusForme service personnel only.**

# SHIATSU & VIBRATION NECK MASSAGER WITH HEAT

## **Heat Button**

Soothing heat for a more relaxing massage.

## **Power Switch**

Turn the unit on by sliding the switch up one notch for Shiatsu Massage, up a second notch for Vibration Massage, and all the way up for Shiatsu and Vibration Massage combined. Slide the switch all the way down to turn the power off.

## **Power Supply Cord**

The Neck Shiatsu and Vibration Massager is powered by a 120 V power cord.

### **Shiatsu Massage**

Shiatsu is a deep kneading circular massage. This unit features a spot Shiatsu massage mechanism that works your neck area.

### **Vibration Massage**

Enjoy invigorating vibration massage or add it to the Shiatsu massage.

### **Convenient and Portable**

Perfect for use at home or in the office.

### **Flex Handles**

Lets you customize the intensity of your massage.

# INSTRUCTIONS FOR USE

To reduce the risk of electrocution, this appliance is equipped with a polarized plug (one blade wider than the other). This plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician to install the proper outlet. Do not change the plug in any way.

Fig. 1

1. Insert the power cord into a 120 volt AC outlet. To activate the Shiatsu massage action, slide switch up one notch. To activate vibration massage, slide switch up a second notch. (Fig. 1)

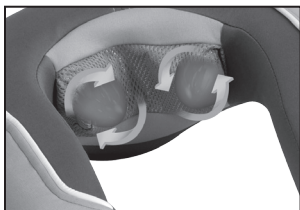
Fig. 2

2. To activate Shiatsu and vibration massage combined, slide switch all the way up. To activate soothing heat with either massage, press the heat button. Massage heads will glow red (Fig. 2). Allow a few moments for the heating action to warm up. When finished slide the switch all the way down and unplug the massager.

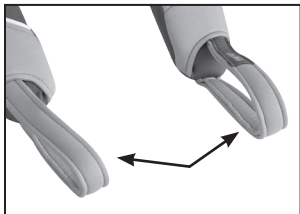
Fig. 3

**NOTE:** For more intense massage, place hands in the comfort grip handles and gently pull forward (see Fig. 3).

# MAINTENANCE

## **To Store**

Place massager in its box or in a safe, dry, cool place. Avoid contact with sharp edges or pointed objects which might cut or puncture the fabric surface. To avoid breakage, do not wrap the power cord around the unit. Do not hang the unit by the hand control cord.

## **To Clean**

Unplug the unit and allow it to cool before cleaning. Spot clean only with a soft, slightly damp cloth. Never allow water or any other liquids to come into contact with the unit.

- Do not immerse in any liquid to clean.
- Never use abrasive cleaners, brushes, gasoline, kerosene, glass / furniture polish or paint thinner to clean.
- Do not attempt to repair the Shiatsu & Vibration Neck Massager. There are no user serviceable parts. For service, call the Consumer Relations telephone number listed in the Warranty Section.

HoMedics Group Canada  
6460 Kennedy Road. Unit C  
Mississauga, Ontario,  
Canada L5T 2X4

Telephone: 1-888-225-7378  
8:30 a.m. to 5:00 p.m.,  
Mon - Fri EST

[www.ObusForme.com](http://www.ObusForme.com)

### **LIMITED WARRANTY (FOR CANADA ONLY)**

ObusForme guarantees this product will be free from defects in workmanship & materials for a period of one (2) years from the date of original purchase. This guarantee applies when the product is used for the purpose intended. In the advent of a defect as described above, please return this product together with the original purchase receipt to ObusForme for refund or replacement. This warranty gives you rights that vary by province/state. This warranty may change.

### **WHAT IS NOT COVERED**

Wear & tear, aging (including foam discoloration, flattening, density, consistency), accidental damages, alterations, mishandling, faulty adjustment, misuse, improper care, power damage, rental use, discontinued items, service by anyone other than ObusForme, shipping damages, neglect, items sold 'as is' or damage due to natural acts are NOT covered.

### **HOW TO OBTAIN WARRANTY SERVICE**

Please obtain a Return Authorization Number & instructions prior to sending your item or it may be denied.

Please contact Customer Service by:

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